measure. don't guess.

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Small and smart: All Medlab products are compact, mobile, robust and can be operated intuitively. Reliability in handheld format or as a desktop model. You are always prepared for safe patient monitoring, regardless of where you are.

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Established in 1993, Medlab is a successful medical device manufacturer in the technology region of Karlsruhe, Germany. High-quality monitoring systems for the measurement and monitoring of vital signs are developed, produced and marketed from our headquarters near Karlsruhe.

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Monitoring-technology of Medlab is the first choice where compact size is essential, but where flexibility is needed at the same time. In hospitals, clinics, sleep labs and in emergency medical services users like the quality of our compact monitoring systems.

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Content

- 02 Medlab
- 03 Monitoring
- 04 Professionals
- 05 Everywhere
- 06 Technology
- 08 Product Overview
- 10 Pulse Oximeters
- 15 Vital Signs Monitors
- 16 Capnographs
- 18 Accessories

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Patient safety is the primary concern in all our developments. Oxygen saturation of the arterial blood, pulse rate, CO₂ concentration of exhaled air, blood pressure and respiration rate are monitored.

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Medlab concentrates in high-end Capnographs, Pulse Oximeters and Vital Signs Monitors. Top-quality for challenging situations: Reliable measurements, clear representation on graphical and numerical displays, progressive acoustical and optical alarms.

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Medlab is one of the world-wide technology-leaders for medical monitoring systems. Every day, we search for new solutions and opportunities. All products are developed and manufactured based on certified quality standards.

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On top of our own product range, Medlab also develops and manufactures modules that are used in medical devices of other world renowned companies.

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Made in the technology region of Karlsruhe, Medlab products are sold world-wide. An extensive distributor network ensures global availability. In some countries, Medlab's medical solutions are market leaders.

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The sales partners of Medlab ensure constant on-site response. Technical information, customer support and product demonstrations are available world-wide. An extensive service-network makes sure that maintenance and repair of Medlab products is executed in a timely manner.

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Intelligent PEARL®-Technology strengthens the signal in pulse oximetry and suppresses artefacts. Advantage: Reliable measurements also during strong patient's movement and low perfusion.

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
The PEARL® algorithm detects SpO₂ and pulse rates reliably, also in difficult situations. The method correlates the measured signals permanently with template curves of different pulse rates and amplifies the correct pulse curve, while attenuating the artefact part of the signal.

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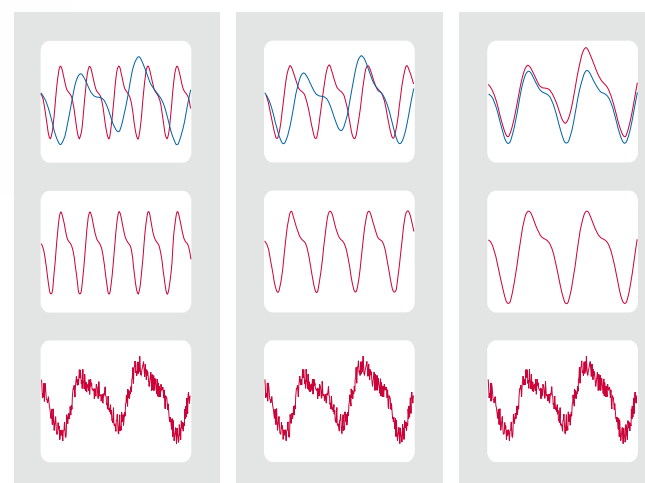PEARL®: Pulse Enhancement Artefact Rejection Logic
PEARL is a registered trademark of Medlab GmbH

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Medlab is certified according to addendum II of the European directive 93/42 EWG, to EN 13485.

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You are on the safe side: All Medlab products, including reusable SpO₂ probes, come with a 24 month warranty. Wear and tear parts, such as battery-packs, have a 6 month warranty.

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Medlab offers you a choice in patient monitoring. You can choose from an extensive range of products that meet all medical and technical requirements.

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* only Capnox
** only PEARL8

Medlab reserves the right to make changes without further notice to any products herein to improve reliability, function or design.

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Lightweight and small-sized Pulse Oximeter for the measurement of pulse rate and oxygen saturation. Trend memory and alarm system (PEARL 8).

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Mobile Pulse Oximeter with batteries. For measurement of pulse rate and oxygen saturation. With alarms and measured data storage of up to 55 hours.

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An optional charger for NANOX10 / PEARL10 expands the mobile Pulse Oximeter to an efficient desktop device.

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Mobile Pulse Oximeter with PEARL® Technology, colour TFT display and data storage of 55 hours. Scope of functions like NANOX10.

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Compact and lightweight Pulse Oximeter for measurement of oxygen saturation, pulse rate and pulse strength. Display of plethysmographic waveform or trended data. Operation by integrated rechargeable battery or mains.

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Lightweight, compact Pulse Oximeter with PEARL® Technology. Measurement of oxygen saturation, pulse rate and pulse strength. Display of plethysmographic waveform or trended data. Operation by integrated rechargeable battery or mains.

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Compact and cost-effective Pulse Oximeter for measurement of oxygen saturation, pulse rate and pulse strength. Display of plethysmogram and numerical values. With trended data storage and full alarms. Operation by integrated rechargeable battery or mains.

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Compact Vital Signs Monitor with non-invasive blood pressure measurement and PEARL® technology supported pulse oximetry. Display of plethysmographic waveform, blood pressure values, pulse, oxygen saturation and trended data. Operation by integrated rechargeable battery or mains.

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Handheld and Desktop

CAPNOS

Capnogram
etCO₂
Respiration frequency

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Light, mobile Sidestream Capnograph with exceptional accuracy.

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Scope of functions like CAPNOS, additionally with PEARL® technology supported pulse oximetry.

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Lightweight, compact sidestream capnograph for measurement of etCO₂ and respiration rate. Display of capnogram or trend. Operation by integrated rechargeable battery or mains.

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An extensive program of accessories completes Medlab's product range. This includes SpO₂ probes, blood pressure cuffs, CO₂ accessories and other accessories such as batteries and protective pouches.

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CO₂ sample lines for intubated and non-intubated patients as well as moisture filters.

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Comprehensive range of high-quality disposable and reusable SpO₂ sensors for devices with Medlab pulse oximetry.

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Holders for attachment to infusion pole or anaesthesia rail for hand and table devices as well as protective pockets and rubber protective covers.

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Blood pressure cuffs in different sizes. Latex-free and antibacterially coated.

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