


EG01010

Three Lead ECG Module

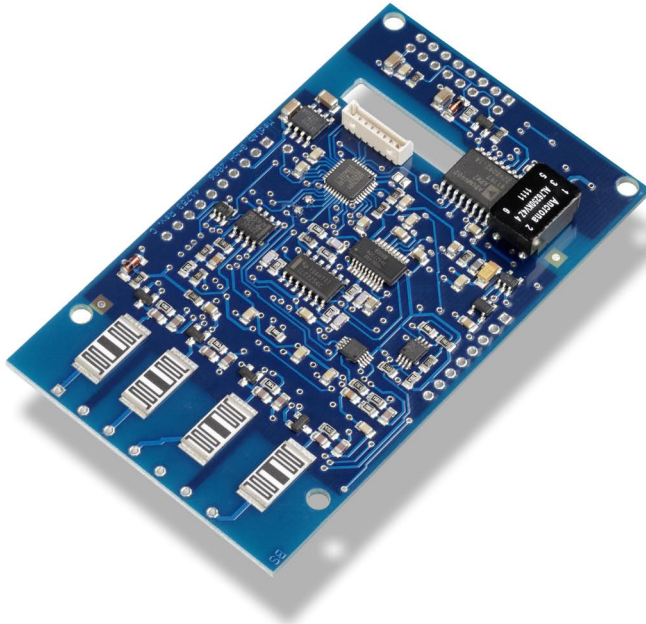

ECG

oem.medlab.eu

- Advanced signal processor technology with digital filtering
- Small size (88 x 53 mm), maximum height 12 mm.
- Monitors Einthoven I or II or III, with three lead cable
- Integrated pacemaker detection and lead-off detection for each electrode
- Measurement of respiration rate over the electrodes (optional)
- Single power supply (5 V DC), < 250 mW
- RS232 and TTL level serial asynchronous interface
- Highly accurate isolated ECG Trigger output
- Integrated patient isolation according to IEC 60601-1
- Integrated defibrillation protection
- Complies with the newest version of IEC 60601-2-27
- Transmission of pulse rate, ECG curve and respiration waveform
- Measurement of pulse rate between 30 – 247 bpm
- Four amplification stages, three curve transmission speeds
- ECG cables and adapters available



Applications:

High-End ECG Monitoring

ECG Triggering for CT scanners, Gamma Cameras etc.

medlab

medizinische Diagnosegeräte GmbH

Helmholtzstrasse 1a

D 76297 Stutensee near Karlsruhe
Germany

Tel.: +49 (0) 7244 74110 0

FAX: +49(0) 7244 74110 28

sales@medlab.eu

oem.medlab.eu