


PEARL200*

Pulse Oximeter OEM Board

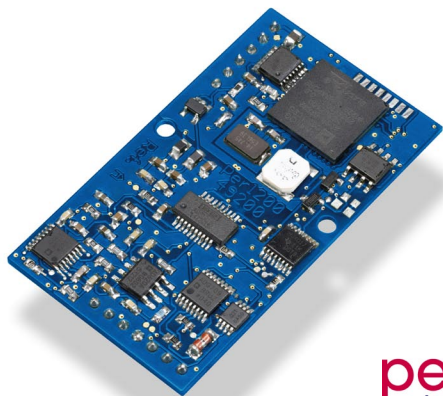

SpO₂

oem.medlab.eu

- Advanced signal processor technology with fourth generation algorithm and artefact suppression, recommended for neonates, infants, and adults
- Very small size (55 x 31.5 mm), maximum height 4 mm
- Also available with pre-installed female or male headers - please specify.
- Position of probe and host connector compatible to PEARL100 board
- Lock-in based signal recovery down to a pulsation of 0.05 % (AC / DC ratio)
- Ultra low noise design for highest performance on centralized patients
- Greatly improved suppression of disturbances caused by ambient light and electromagnetic interference
- Easily integrated into handheld and desktop devices
- Single power supply (3.3 V DC), < 100 mW for average finger thickness (20 – 40mA, depending on LED brightness, adapts to finger properties)
- TTL level serial asynchronous interface (signal levels 0 V and 3.3 V)
- Transmission of SpO₂, pulse rate, plethysmographic waveform, perfusion index, "quality"-signal, and status information
- Measurement of pulse rate between 30 and 250 bpm
- Measurement of oxygen saturation between 0 and 100 %
- Complete line of probes available
- Complies with applicable clauses of ISO 9919 and ISO 80601-2-61



medlab

medizinische Diagnosegeräte GmbH

Helmholtzstrasse 1a

D 76297 Stutensee near Karlsruhe

Germany

Tel.: +49 (0) 7244 74110 0

FAX: +49(0) 7244 74110 28

sales@medlab.eu

oem.medlab.eu

* Pulse Enhancement Artifact Rejection Logic
PEARL is a registered trademark of Medlab GmbH